

VISALYS® BULK FLOW


KETTENBACHDENTAL
Simply intelligent



**STRIKINGLY
ESTHETIC**





FLOWABLE ONE-SHADE BULK COMPOSITE

VISALYS® BULK FLOW

VISALYS® BULK FLOW, our light curing, flowable nano-hybrid filling composite for filling cavities up to 4 mm depth. It combines efficiency, stability and esthetics.

ONE-SHADE
SYSTEM

Fascinating shade adaption with just one shade thanks to the One-Shade system. A real chameleon effect with natural fluorescence for restorations close to nature.

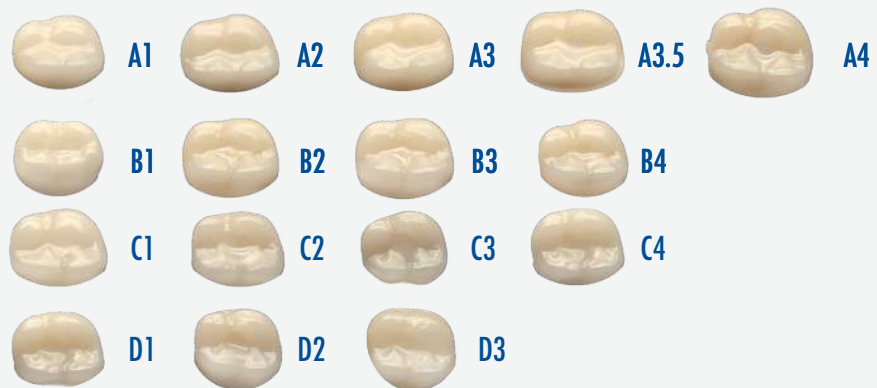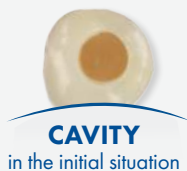

Exemplary fillings using VISALYS® BULK FLOW on Vita Shades



Excellent polishability for high gloss results

USER FEEDBACK

VISALYS® BULK FLOW convinces with a pleasant handling, a good color adaptation and an excellent gloss.

Dr. Maximilian Dobbertin

CONVINCING CONSISTENCY

VISALYS® BULK FLOW behaves like a flow but is stable and modelable like a fill.

It provides high stability and at the same time high radiopacity. Low shrinkage stress ensures a safe adaption to the cavity walls.

No oozing thanks to an intelligent syringe design.

Precise fillings without wasting material thanks to the build in movable seal.




CONSISTENCY
OF A
FLOW
WITH THE
MODELABILITY
OF A

FILL





STRIKINGLY ESTHETIC

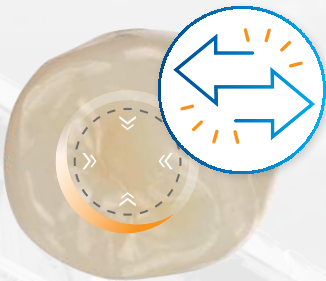

The **real chameleon effect** of our **One-Shade System** and a **natural fluorescence** enable vivid, natural-looking restorations.



High covering capacity or natural aesthetic results already from a depth of cure of 3 mm.*



VISALYS® BULK FLOW



OMNICHROMA FLOW BULK

*Exemplary 3 mm discs fabricated from the respective filling materials, shade „Universal“



Efficient working **without increments** or **additional capping layer**



**FREE FROM
BIS-GMA**




VERY DENSE MARGINS

Excellent adaption to the cavity walls reduces the risk of marginal gaps and marginal discoloration for long lasting restorations.


Marginal quality
on enamel

Overview chewing simulator with Marburg chewing simulator (TCML)*

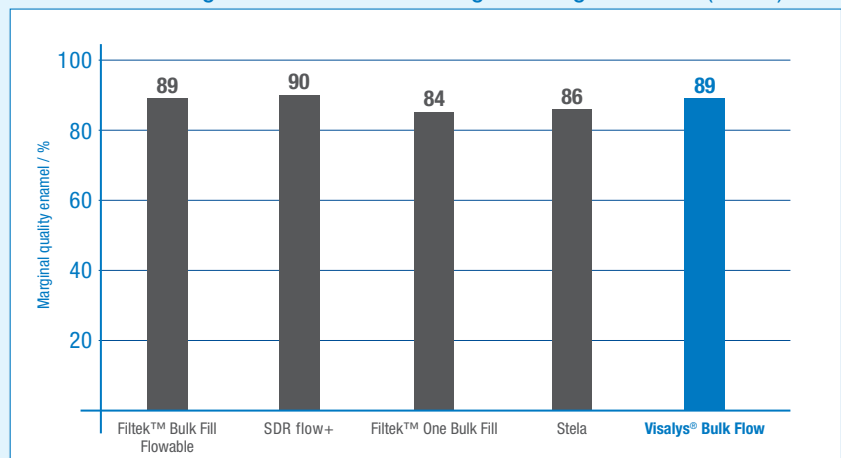

Marginal quality
on dentin


Overview chewing simulator with Marburg chewing simulator (TCML)*

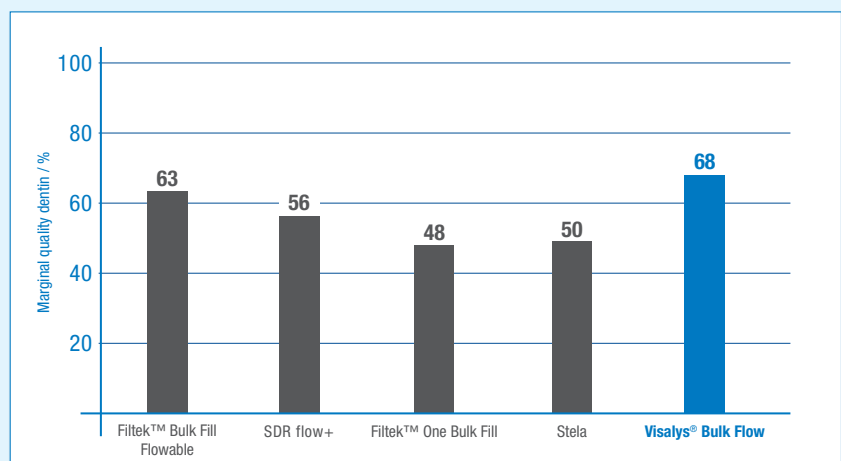

* Messungen Prof. Roland Frankenberger, Universität Marburg, 2024.
Thermo-Mechanische Wechselbelastung: 1.000.000 Zyklen à 50 N, 25.000 thermal cycles (5/55 °C), 264 teeth (n=8), MOD class II.

THE RANGE

VISALYS® BULK FLOW

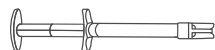

1 x 2 g Syringe
5 x Application tips

REF 15064

VISALYS® BULK FLOW



15 x 0,25 g Caps

REF 15066

2 x 0,25 g Caps

REF 15067



Made in
Germany

Kettenbach GmbH & Co. KG
Im Heerfeld 7
35713 Eschenburg · Germany
Phone +49 (0)2774 7050
info@kettenbach.com
www.kettenbach-dental.com



DEN1410933_1126

KETTENBACHDENTAL
Simply intelligent