




SO ADHERENT IS **MUCOPREN®** SOFT

SOFT RELINING

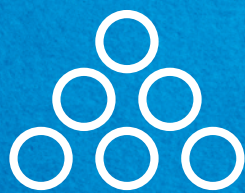


KETTENBACHDENTAL

Simply intelligent

MUCOPREN® SOFT

Prosthetic relining

Mucopren® Soft is a soft long-term relining material based on A silicone for **direct and indirect application**.

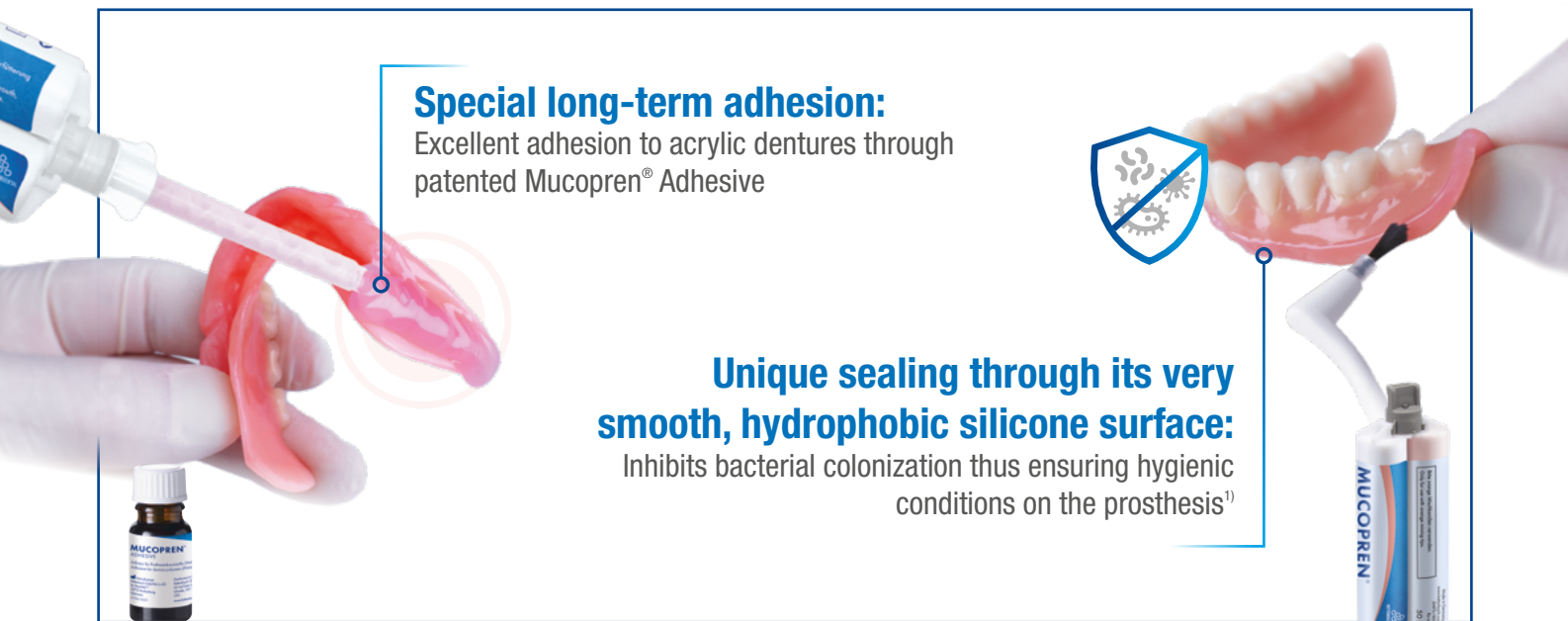

Special long-term adhesion:
Excellent adhesion to acrylic dentures through patented Mucopren® Adhesive

Unique sealing through its very smooth, hydrophobic silicone surface:
Inhibits bacterial colonization thus ensuring hygienic conditions on the prosthesis¹⁾



Step by Step

Simple treatment & fast processing

Faster setting times and shorter waiting times ensure a fast and smooth process in the practice and laboratory.

Supports the healing process of implants

Also suitable for use in implantology for trouble-free osseointegration.



High wearing comfort

The definitive resilient reline is completely odorless and tasteless and makes wearing the prosthesis as pleasant as possible.



HOLDS TO ITS PROMISES



FURTHER INDICATIONS



- For the creation of **permanently soft relinings**
- For relieving pressure on the surgical site after insertion of **implants**
- For balancing mucous membrane areas with different resilience
- For relieving sharp bony ridges
- For relining **interim prostheses**
- As **mucous membrane conditioner** in case of prosthetic stomatitis
- For creating and relining **obturators**

Benefits of the optimized adhesive²⁾

» ONLY A **SINGLE APPLICATION** NECESSARY



3 MIN

CHAIRSIDE

The processing time is reduced to only **3 minutes** following mouth removal. No water bath necessary.



3 MIN

LABSIDE

The processing time is reduced to only **3 minutes** in the water bath or pressure pot.

¹⁾ Result from in-vitro studies in which it was shown that the growth of *Candida albicans* was prevented thanks to the smooth silicone surface

²⁾ Refers to the optimized product version of Mucopren® Adhesive, first introduced in 2023.

MUCOPREN® SOFT – THE RANGE



Mucopren® Soft

Mucopren® Soft Basic Set	REF 28105	50 mL Mucopren® Soft, 50 mL Mucopren® silicone sealant, 10 mL Mucopren® Adhesive, 6 mixing tips white, 16 mixing tips orange, 1 brush holder, 20 single-use brushes, 1 steel bur, accessories
Mucopren® Soft Normal pack	REF 15687	2 x 50 mL Mucopren® Soft, 6 mixing tips white
Mucopren® silicone sealant	REF 15686	50 mL Mucopren® silicone sealant, 8 mixing tips orange
Mucopren® Adhesive	REF 14203	10 mL Mucopren® Adhesive



Technical data

Product name	Mixing ratio	Working time at 23 °C / 74 °F ≤	Intraoral setting time at 35 °C / 95 °F ≥	Water bath	End of setting ≥	Shore hardness (approx.)
Mucopren® Soft – direct relining	1:1	2:15 min.	3:15 min.	Not applicable	5:30 min.*	A 28
Mucopren® Soft – indirect relining	1:1	2:15 min.	Not applicable	3:00 min.	> 3 min.	A 28
Mucopren® silicone sealant	1:1	1:30 min.	Not applicable	Not applicable	6:30 min. at 23 °C / 74 °F	A 33

*Total setting time (removal time from mouth) from start of mixing

NOTES

Kettenbach GmbH & Co. KG
Im Heerfeld 7
35713 Eschenburg · Germany
Telefon +49 (0)2774 7050
info@kettenbach.com
www.kettenbach-dental.com



1410871_4924

KETTENBACHDENTAL
Simply intelligent